

Owner's Manual

Infrared Ear/Forehead Thermometer

Model DET-2126



Importer:



JOYTECH Healthcare Co., Ltd.
No.365, Wuzhou Road,311100 Hangzhou, Zhejiang Province,
PEOPLE'S REPUBLIC OF CHINA

Telephone: +86-571-81957767
Fax: +86-571-81957750



The product is in compliance with the requirements of
MDR(EU) 2017/745, "0123" is the identification number
of notify body



European Authorized Representative:
Shanghai International Holding Corp. GmbH (Europe)
Eiffelstrasse 80, 20537 Hamburg, Germany

Document No.: JDET-3004-040
Article No.: 003
Version:A
Date of Issue: 2025-4-21

Contents

Introduction	02
Product Description	03
Lcd Display Introduction	04
Basic Functions	05
Body Temperature	06
Illustration For Use	07
Temperature Taking Hints	14
Memory Mode	16
Care And Cleaning	17
Battery Replacement	19
Specifications	20
Troubleshooting	21
Calibration	23
Service	24
Warranty	25
Electromagnetic Compatibility Information	26

Introduction

The ear/forehead thermometer is intended for measuring forehead temperature and ear temperature. It allows you to quickly and easily measure body temperature. It converts the measured heat into a temperature reading displayed on the LCD.

Intended use: The infrared ear/forehead thermometer is intended for the intermittent measurement of human body temperature. The device can be reused by people of all ages for home use and clinical use.

Indications for Use: The infrared ear/forehead thermometer is used to measure body temperature.

Intended user: People of all ages. There is no gender or age limitation for using infrared ear/forehead thermometers.

Claims: No specific claim except the achievement of the intended use.
Complications: None.




Please read all instructions carefully and thoroughly before using this product.

Warning:

1. It's not meant to replace a visit to the doctor. Please consult with doctor if you have health concerns.
2. There is no gender or age limitation for using infrared ear/forehead thermometer.
3. Do not allow children to take their temperatures unsupervised, some parts are small enough to be swallowed.
4. Never immerse this device in water or other liquids (not waterproof).
5. Do not modify this equipment without authorization of manufacturer.
6. Do not expose the thermometer to temperature extremes (below -25°C/-13°F or over 55°C/131°F) nor excessive humidity (>95%RH).
7. Keep the battery away from children.
8. Remove battery from the device when not in operation for a long time.
9. Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.
10. ME equipment should not be cleaned while in use.
11. The probe of the ME equipment shall not be serviced or maintained while in use with a patient.


Product Description

1. Probe Cover
2. Probe
3. Test Button/Power Button
4. Memory Button/Age Button
5. Setting Button/Sound Button
6. Battery Cover

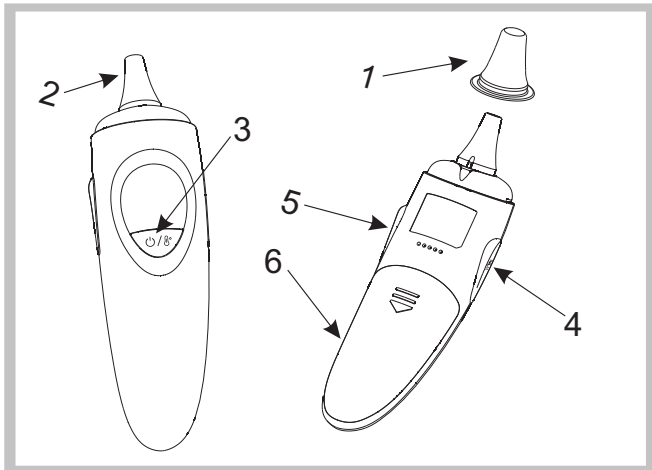

3

LCD Display Introduction

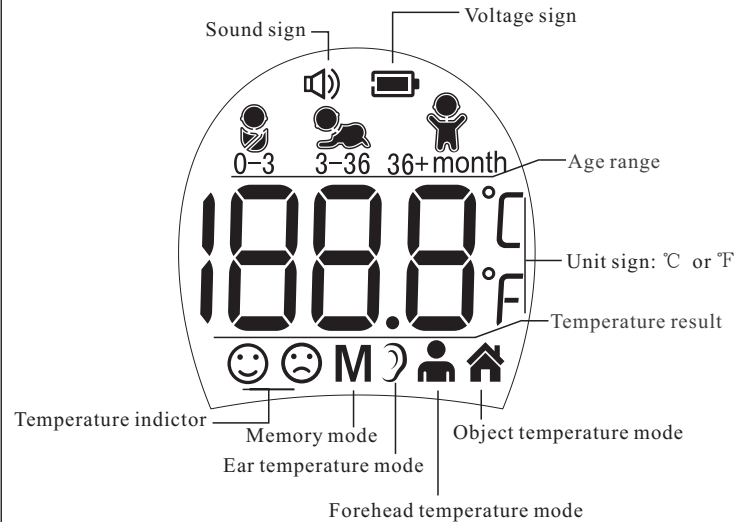

4

Basic Functions

Ear/Forehead Temperature Mode	The thermometer has been designed for practical use. It's not meant to replace a visit to the doctor. Please also remember to compare the measurement result to your regular body temperature. → Please see the Illustration For Use section to learn how to measure the body temperature.
Object Temperature Mode	The object mode shows the actual, unadjusted surface temperatures, which is different from the body temperature. It can help you to monitor if the object temperature is suitable for the baby or patient, for example the baby's milk. Measure time of object mode: Approx. one second. Measuring range of object mode: 0°C~100°C (32°F~212°F). Laboratory accuracy of object mode: ±4% or ±2°C (4°F) whichever is greater. → Please see the Illustration For Use section to learn how to measure the object temperature.
Memory Mode	There are 30 sets memories for ear, forehead and object measurements. Each memory also records the measurement mode and age range icon.
°C/°F Switch	Please see Illustration For Use section to learn how to change between Celsius and Fahrenheit.
Beep Hint	When the testing is done, it will sound one short beep.
Vibration	After the beep is turned off, the vibration will be automatically turned on. The thermometer will vibrate at the end of the measurement.
Age range	Please see Illustration For Use section to learn how to select age range among 0-3 months, 3-36 months and 36 months-adult.

5

Body Temperature

The temperature of a healthy person is affected by various factors: the person's individual metabolism, their age (body temperature is higher in babies and toddlers and decreases with age. Greater temperature fluctuations occur faster and more often in children, e.g. due to growth spurts), their clothing, the ambient temperature, the time of day (body temperature is lower in the morning and increases throughout the day towards evening), the preceding physical and, to a lesser extent, mental activity.

It is suggested that users develop the habit of daily measurement, and use this product to establish a private normal temperature range for themselves and their families as a reference for measurement.

6

Illustration For Use

► Before measurement:

- Check before each use that the lens is intact. If it is damaged, contact your retailer or the service address.
- Bear in mind that the thermometer needs to have been in the room in which the measurement is taken for at least 30 minutes before use.
- Users and the thermometer should be in the same ambient temperature.
- Before each measurement, make sure that the device is in the appropriate mode for the measurement that you wish to take.

► How to switch between ear/forehead and object mode:

- When the device is on, press and hold the *Setting Button/Sound Button* to switch between ear mode/forehead mode and object mode. When the forehead cap is installed, the thermometer will automatically switch to forehead mode and when not installed, the thermometer will automatically switch to ear mode.

► To measure forehead temperature:

Note:

- Oils or cosmetics on the forehead may give a lower temperature reading than the actual one. Remove dirt from the forehead before taking a measurement. Wait at least 10 minutes after washing the forehead area before taking a reading.

10

Illustration For Use

- Do not use the thermometer on a perspiring or sweating forehead, as this may affect the reading.
 - Remove hat and hair and wait 10 minutes before taking a reading.
1. Install the forehead cap.
 2. Press the *Test Button/Power Button* to switch on the thermometer and the display is activated to show all segments.
 3. After self-checking, Figure 1 appears on the display screen with beep, so you can start a new measurement.
Note: You can press Memory Button/Age Button to change age range among 0-3 months, 3-36 months and 36 months-adult.
 4. Place the probe with the forehead cap against the forehead (See Figure 2).
 5. Press the *Test Button* to start measurement.
 6. The measured value appears on the display and the end of the measuring is signaled by beep. If the beep is turned off, the end of the measurement is signaled by vibration.
 7. The thermometer will shut off automatically after 30 seconds of inactivity. To prolong battery life, press and hold the *Test Button* to turn the unit off.

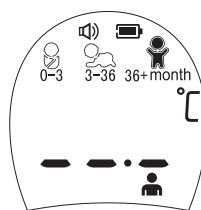


Figure 1



Figure 2

8

Illustration For Use

► To measure ear temperature:

- Some people produce different readings in their left and right ear. In order to record temperature changes, always measure a person's temperature in the same ear.
 - As ear wax can affect the measurement, you should clean the ear before measuring if necessary.
1. Remove the forehead cap.
 2. Press the *Test Button/Power Button* to switch on the thermometer and the display is activated to show all segments.
 3. After self-checking Figure 3 appears on the display screen with beep, so you can start a new measurement.
 4. Make sure that the sensor tip and the ear canal are clean. As the ear canal is slightly curved, you have to pull the ear slightly up and backwards before inserting the sensor tip. (See Figure 4)
 - When using the thermometer on infants under age 1, pull the ear up making sure the sensor faces the eardrum (See Figure 5).
 - When using the thermometer on individuals over the age of 1, pull the ear back making sure the sensor faces the eardrum. (See Figure 6).
5. Press the *Test Button/Power Button* to start the measurement.

9

Illustration For Use

6. The measured value appears on the display and the end of the measuring is signaled by beep.
Note: A waiting period of 20 seconds between testing is recommended to avoid excessive cooling of the skin.
7. The thermometer will shut off automatically after 30 seconds of inactivity. To prolong battery life, press and hold the *Test Button/Power Button* to turn the unit off.

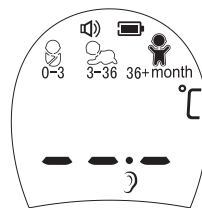


Figure 3

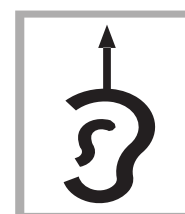


Figure 5



Figure 6

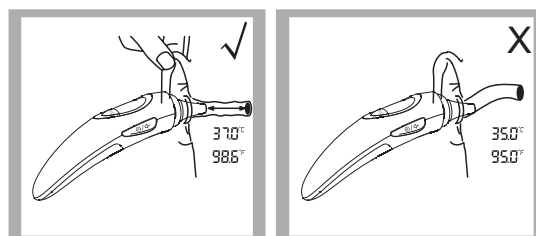


Figure 4

Illustration For Use

► To measure object temperature:

1. Press the *Test Button/Power Button* to switch on the thermometer.
2. Press and hold *Setting button/Memory Button* to switch to object temperature mode. Figure 7 appears on the display screen with beep, so you can start a new measurement.
3. Aim the thermometer at the center of the object you want to measure with a distance of 1 to 2 cm.
4. Press the *Test Button/Power Button* to start measurement.
5. The measured value appears on the display and the end of the measuring is signaled by beep. If the beep is turned off, the end of the measurement is signaled by vibration.
6. To exit object temperature mode, press and hold *Setting Button/Memory Button* to switch back to forehead/ear mode.
7. The thermometer will shut off automatically after 30 seconds of inactivity. To prolong battery life, press and hold the *Test Button/Power Button* to turn the unit off.

Remarks:

This mode shows the actual, unadjusted surface temperatures, which is different from the body temperature.

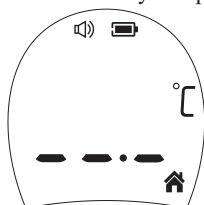


Figure 7

11

Illustration For Use

► How to turn on or off sound:

You can press *Setting Button/Memory Button* to turn on or off sound with thermometer on.

► How to set unit between °C and °F:


With the thermometer off, press and hold *Setting Button/Memory Button* to enter into setting mode. Press *Test Button/Power Button* to select the unit you want. After the unit is set, press *Setting Button/Memory Button* to exit the setting mode (See Figure 8).



Figure 8

► How to select the age range:

With the thermometer on, press the *Memory Button/Age Button* to select the age range among 0-3 months, 3-36 months and 36 months-adult (See Figure 9).

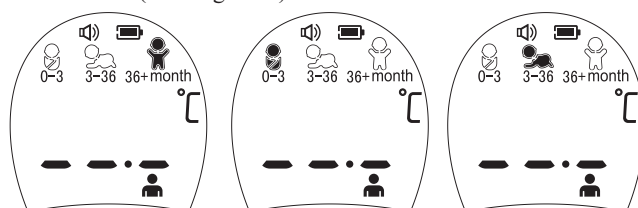

► After measurement: Figure 9

1. Power off: Device will automatically shut off if left idle for more than 30 seconds to extend battery life.
2. Clean the probe after each use to ensure an accurate reading and avoid cross contamination.
(See the section of Care and Cleaning for details.)

12

Illustration For Use

► BackLight:

In Ear/Forehead mode:

The display will be lighted GREEN for 3 seconds when a measurement is completed.

Age range	GREEN 😊	YELLOW 😊	RED 😊
0-3 Months	34.0°C-37.4°C (93.2°F-99.4°F)		37.5°C-43.0°C (99.5°F-109.4°F)
3-36 Months	34.0°C-37.6°C (93.2°F-99.6°F)	37.7°C-38.5°C (99.7°F-101.3°F)	38.6°C-43.0°C (101.4°F-109.4°F)
36 Months-adult	34.0°C-37.7°C (93.2°F-99.9°F)	37.8°C-39.4°C (100.0°F-103.0°F)	39.5°C-43.0°C (103.1°F-109.4°F)

In Object mode:

The display will only be lighted Green for 3 seconds when the unit is ready for measurement and a measurement is completed.

13

Temperature Taking Hints

To ensure that the reading always reflects the body temperature accurately, you need to take account of the following factors which may affect an accurate reading.

1. It is important to know each individual's normal temperature when they are well. This is the only way to accurately diagnose a fever. To determine normal temperature, take multiple readings when healthy. Re-measure with a standard digital thermometer for confirmation.
2. The probe window of the thermometer is the most delicate part of the device. Do not touch the probe window. The accuracy of the reading may be affected if the probe window is damaged or dirty.
3. Users should not drink, eat, or be physically active such as bathing, showering, shampooing and hair drying before/while taking the measurement. In these cases, wait 20 minutes prior to taking a temperature.
4. Holding a hand on the forehead for any length of time will affect the temperature reading.

14

Temperature Taking Hints

- Do not take temperature over scar tissue, open sores or abrasions.
- Don't take a measurement while or immediately after nursing a baby.
- Do not use this thermometer outdoors.
- Do not take temperatures with this thermometer near places that are very hot, such as fireplaces and stoves.
- If the thermometer is stored in a significantly different environment than testing location, place it in the testing location for approximately 30 minutes prior to use.
- The patient is an intended operator.

15

Memory Mode

- With the thermometer off, press the *Memory Button/Age Button* to enter memory mode, the letter M will appear in the upper left corner of the display. (See Figure 10)
- The thermometer will automatically memorize the last 30 temperature readings. Each memory also records the measurement mode and age range icon. Each time the *Memory Button/Age Button* is pressed, the screen displays past readings that correspond with a number 1-30. The number 1 reflects the most recent reading, while the number 30 reveals the oldest reading stored in memory. (See Figure 11)
- In the memory mode, the mode and age range icon always exists. The user can press the *Test Button/Power Button* to take new measurements.

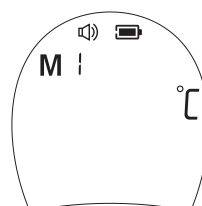


Figure 10

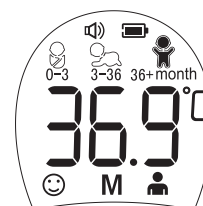


Figure 11

16

Care And Cleaning

- The probe window must be kept clean, dry, and undamaged at all times to ensure accurate readings. The accuracy of temperature readings can be affected by damage to the probe window, or the presence of dirt and ear wax on the probe window.
- Fingerprints, earwax, dust and other soiling compounds reduce transparency of the window and result in lower temperature readings.
- The probe window is the most delicate part of the thermometer. To safely clean the window, gently wipe its surface with a cotton swab slightly moistened with isopropyl alcohol and immediately wipe dry with a clean cotton swab. After cleaning, allow at least 5 minutes drying time before taking temperatures.
Note: Do not use any chemical other than isopropyl alcohol to clean the probe window.
- Use a soft, dry cloth to clean the thermometer display and exterior.
- The thermometer is not waterproof. Do not submerge the unit in water when cleaning.
- Store the thermometer in a dry location, free from dust and contamination and away from direct sunlight.

17

Care And Cleaning

- Periodic cleaning of the device following use to prevent patient cross infection.
-Use a soft cloth slightly moistened with a 75% isopropyl alcohol solution to wipe the thermometer and probe. Do not use abrasive cleaners.
- Ensure that children do not use the instrument unsupervised; some parts are small enough to be swallowed.
- Do not remove or modify the equipment without permission.
- Strong electromagnetic fields may interfere with the proper operation of the thermometer. The device needs special pre-cautions regarding EMC according to the EMC information.
- It is not intended for use in the oxygen rich environment and presence of flammable anesthetic mixture with air, oxygen or nitrous oxide.
- Put the thermometer back to the original packaging after using.

18

Battery Replacement

1. Replace battery when "  " appears in the upper right corner of LCD display. (See Figure 12)
2. Slide battery cover down as shown in Figure 13.
3. Remove battery and install 2 new AAA alkaline batteries as shown in Figure 14.
4. Slide battery cover back on.

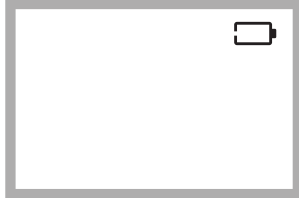


Figure 12

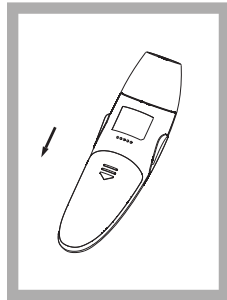


Figure 13

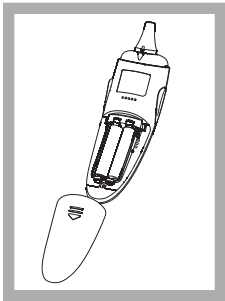


Figure 14

19

Specifications

Measuring range	Ear/Forehead mode: 34.0°C~43.0°C(93.2°F~109.4°F)
Measuring site	Ear canal(Ear Mode)、Forehead(Forehead Mode)
Reference body site	Oral (This thermometer converts the ear temperature to display its "oral equivalent.")
Operation mode	Ear/Forehead mode(Adjust mode)
Laboratory accuracy	Ear/Forehead mode: ±0.2°C (0.4°F) during 35.5°C~42.0°C (95.9°F~107.6°F) at 15°C~35°C (59.0°F~95.0°F) operating temperature range ±0.3°C (0.5°F) for other measuring and operating temperature range
Display resolution	0.1°C or 0.1°F
Measure time	Ear/Forehead mode: Approx. one second
Operating temperature range:	10°C~40°C(50°F~104°F), 15%~85%RH, non-condensing Atmospheric Pressure : 70kPa ~ 106kPa
Storage and transport temperature range	-25°C ~ 55°C (-13°F ~ 131°F), 15%~95%RH, non-condensing Atmospheric Pressure : 70kPa ~ 106kPa
Clinical accuracy	0-1 year: Ear temperature: Clinical bias: 0.05°C (0.09°F) ; Clinical repeatability: 0.21°C (0.38°F) ; Limits of agreement: 0.76°C (1.37°F) Forehead temperature: Clinical bias: 0.08°C (0.14°F) ; Clinical repeatability: 0.33°C (0.59°F) ; Limits of agreement: 0.73°C (1.31°F) 1-5 years: Ear temperature: Clinical bias: 0.11°C (0.20°F) ; Clinical repeatability: 0.34°C (0.61°F) ; Limits of agreement: 1.36°C (2.45°F) Forehead temperature: Clinical bias: 0.24°C (0.43°F) ; Clinical repeatability: 0.40°C (0.72°F) ; Limits of agreement: 1.49°C (2.68°F) over 5 years: Ear temperature: Clinical bias: 0.21°C (0.38°F) ; Clinical repeatability: 0.21°C (0.38°F) ; Limits of agreement: 0.99°C (1.78°F) Forehead temperature: Clinical bias: 0.22°C (0.40°F) ; Clinical repeatability: 0.26°C (0.47°F) ; Limits of agreement: 1.10°C (1.98°F)
Shock	withstands drop of 3 feet
Dimension	155*43*30.5mm
Weight	Approx.93grams(with batteries)
Battery	DC3V(2×AAA battery)
Battery life	Approx. 6000 readings
Expected service life/ Shelf life	Three years
Ingress protecting rating	IP22
Accessories	Battery
Contraindications	None
Applied part	Probe

20

Troubleshooting

Error message	Problem	Solution
---	Measurement before thermometer is ready	Take a measurement until thermometer is ready(with beep).
Er-2	The ambient temperature is not within the range between 10°C and 40°C (50°F~104°F).	Place the thermometer in a room for at least 30 minutes at room temperature between 10°C and 40°C (50°F~104°F)
Er-3	The thermometer is placed incorrectly or unsteady.	Read Illustration For Use thoroughly and take a new temperature measurement.
Er-4	The thermometer showing a rapid ambient temperature change.	Allow the thermometer to rest in a room for at least 30 minutes at room temperature: between 10°C and 40°C (50°F~104°F)
Er-5	The thermometer is not functioning properly.	Unload the battery, wait for 1 minute and repower it. If the message reappears, contact the retailer for service.

21

Troubleshooting

Error message	Problem	Solution
H _i	In Ear/Forehead mode: Temperature taken is higher than 43.0 °C (109.4°F). In Object mode: Temperature taken is higher than 100 °C (212°F).	Read Temperature Taking Hints Thoroughly, then take a new temperature measurement.
Lo	In Ear/Forehead mode: Temperature taken is lower than 34.0 °C (93.2°F). In Object mode: Temperature taken is lower than 0°C(32°F).	Read Temperature Taking Hints thoroughly , then make sure the lens filter are clean, then take a new temperature measurement.
	The thermometer works properly.	Use the thermometer normally
	When battery outline flashes, it indicates that the power is low, but you can continue to measure.	The thermometer will take a proper measurement but batteries must be replaced soon.
	The thermometer could not work due to low battery.	Replace two new alkaline batteries size AAA.

22

Calibration

The thermometer is initially calibrated at the time of manufacture. If the thermometer is used according to the use instruction, periodic readjustment is not required. However, We recommends checking calibration every year or whenever clinical accuracy of the thermometer is in question. Please send the complete device to the dealers or manufacturer.

The above recommendations do not supersede the legal requirements. The user must always comply with legal requirements for the control of the measurement, functionality, and accuracy of the device which are required by the scope of relevant laws, directives or ordinances where the device is used.

A clinical summary and procedures for checking calibration are available upon request. (Turn on the thermometer and press the test button long time until entering into calibrate mode, software version will be displayed.)



TYPE BF APPLIED PART



Unique device identifier



Disposal of this product and used batteries should be carried out in accordance with the national regulations for the disposal of electronic products



Refer to instruction manual/booklet



Batch Code



Storage and Transportation Temperature Limit:
-13°F ~ 131°F (-25°C ~ 55°C)



Atmospheric pressure limitation



Medical Device



Date of manufacture



Storage and Transportation Humidity limitation: 15%~95%RH

23

Service

The thermometer has a limited one year warranty. Do not attempt to disassemble or repair the thermometer by yourself. Should service be required during or after the warranty period you must contact the manufacturer. Repackage the thermometer carefully in its original packaging or securely pack to avoid damage during shipping. Include the original sales slip indicating the date of purchase, a note describing the problem, and your return address. Send the thermometer prepaid and insured.

The lay operator or lay responsible organization should contact the manufacturer or the manufacturer's representative:

- for assistance, if needed, in setting up, using or maintaining the thermometer; or
- to report unexpected operation or events.



Manufacturer



General symbol for recovery/recyclable



The product conforms to the requirements of the EC Directive MDR (EU) 2017/745 on medical devices



European Union Authorized representative



The first num.2: Protected against solid foreign objects of 12,5 mm Ø and greater.
The second num.2: Protection against vertically falling water drops when ENCLOSURE tilted up to 15°



Caution

24

Warranty

This appliance conforms to the following standards:

ISO 80601-2-56 Medical electrical equipment —Part 2-56: Particular requirements for basic safety and essential performance of clinical thermometers for body temperature measurement,
IEC 60601-1-11 Medical electrical equipment —Part 1-11: General requirements for basic safety and essential performance –Collateral Standard: Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment and complies with the requirements of IEC 60601-1-2(EMC) ,
IEC 60601-1(Safety) standards.

Thermometer is warranted by manufacture to be free from defects in material and workmanship under normal use and service for a period of one year from the date of delivery to the first user who purchases the instrument. This warranty does not cover batteries, damage to the probe window, or damage to the instrument caused by misuse, negligence or accident, and extends to only to the first purchaser of the product.

25

Electromagnetic Compatibility Information

The device satisfies the EMC requirements of the international standard IEC 60601-1-2. The requirements are satisfied under the conditions described in the table below. The device is an electrical medical product and is subject to special precautionary measures with regard to EMC which must be published in the instructions for use. Portable and mobile HF communications equipment can affect the device. Use of the unit in conjunction with non-approved accessories can affect the device negatively and alter the electromagnetic compatibility. The device should not be used directly adjacent to or between other electrical equipment.

26

Electromagnetic Compatibility Information

Table 1

Guidance and manufacturer’s declaration – electromagnetic emission		
The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment – guidance
Conducted emission CISPR 11	Not applicable	The device uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
Radiated emission CISPR 11	Group 1 Class B	
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations / flicker emissions IEC 61000-3-3	Not applicable	

Electromagnetic Compatibility Information

Table 2

Guidance and manufacturer’s declaration – electromagnetic immunity			
The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	± 8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrostatic transient / burst IEC 61000-4-4	± 2 kV for power supply lines 100 kHz repetition frequency ± 1 kV for input/output lines	N/A	N/A
Surge IEC 61000-4-5	± 0.5 kV, ± 1 kV differential mode line-line	N/A	N/A
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT (100 % dip in UT) for 0.5 cycle at 0°, 45°, 90°, 135°, 180°, 225°, 270°, and 315° 0 % UT (100 % dip in UT) for 1 cycle at 0° 70 % UT (30 % dip in UT) for 25/30 cycles at 0° 0 % UT (100 % dip in UT) for 250/300 cycle at 0°	N/A	N/A
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m, 50/60Hz	30 A/m, 50/60Hz	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Proximity magnetic fields	IEC 61000-4-39	See Table 3	See Table 3
NOTE: UT is the a. c. mains voltage prior to application of the test level.			

Electromagnetic Compatibility Information

Table 3

Test specifications for ENCLOSURE PORT IMMUNITY to proximity magnetic fields		
Test frequency	Modulation	IMMUNITY TEST LEVEL (A/m)
30 kHz a)	CW	8
134,2 kHz	Pulse modulation b) 2,1 kHz	65 c)
13,56 MHz	Pulse modulation b) 50 kHz	7,5 c)
a) This test is applicable only to ME EQUIPMENT and ME SYSTEMS intended for use in the HOME HEALTHCARE ENVIRONMENT. b) The carrier shall be modulated using a 50% duty cycle square wave signal. c) r.m.s., before modulation is applied.		

Table 4

Guidance and manufacturer’s declaration – electromagnetic immunity		
The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.		
Immunity test	IEC 60601 test level	Compliance level
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz 6 Vrms 150 kHz to 80 MHz outside ISM bandsa	N/A
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.7 GHz	10 V/m

Electromagnetic Compatibility Information

Table 4 continued

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.
a The ISM(industrial, scientific and medical) bands between 0,15 MHz and 80 MHz are 6,765 MHzto6,795 MHz; 13,553 MHz to 13,567 MHz;26,957 MHz to 27,283 MHz; and 40,66 MHz to 40,70 MHz. The amateur radio bands between 0,15 MHz and 80 MHz are 1,8 MHz to 2,0 MHz,3,5 MHz to 4,0 MHz,5,3 MHz to 5,4 MHz,7 MHz to 7,3 MHz,10,1 MHz to 10,15 MHz,14 MHz to 14,2 MHz,18,07 MHz to 18,17 MHz,21,0MHz to 21,4MHz,24,89 MHz to 24,99 MHz,28,0 MHz to 29,7 MHz and 50,0 MHz to 54,0 MHZ.
b The compliance levels in the ISM frequency bands between 150 kHz and 80 MHz and in the frequency range 80 MHz to 2,7 GHz are intended to decrease the likelihood that mobile/portable communications equipment could cause interference if it is inadvertently brought into patient areas. For this reason, an additional factor of 10/3 has been incorporated into the formulae used in calculating the recommended separation distance for transmitters in these frequency ranges.
c Field strengths from fixed transmitters, such as base stations for radio(cellular/ cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the device is used exceeds the applicable RF compliance level above, the device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the device.
d Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

Electromagnetic Compatibility Information

Table 5

Recommended separation distances between RF wireless communications equipment				
The device is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the device can help prevent electromagnetic interference by maintaining a minimum distance between RF wireless communications equipment and the device as recommended below, according to the maximum output power of the communications equipment.				
Frequency MHz	Maximum Power W	Distance	IEC 60601 Test Level	Compliance Level
385	1.8	0.3	27	27
450	2	0.3	28	28
710	0.2	0.3	9	9
745				
780				
810	2	0.3	28	28
870				
930				
1720	2	0.3	28	28
1845				
1970				
2450	2	0.3	28	28
5240	0.2	0.3	9	9
5500				
5785				

Note 1: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Electromagnetic Compatibility Information

WARNINGS!

- This device should not be used in the vicinity or on the top of other electronic equipment such as cell phone, transceiver or radio control products. If you have to do so, the device should be observed to verify normal operation.
- The use of accessories and power cord other than those specified, with the exception of cables sold by the manufacturer of the equipment or system as replacement parts for internal components, may result in increased emissions or decreased immunity of the equipment or system.
- Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation.
- Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.
- Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the device, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.
- Portable and mobile RF communications can affect the device. The device needs special pre-cautions regarding EMC according to the EMC information provided in the accompany documents.
- Do not use the devices in the MR environment.
- The Operator should not use the system and should inform the customer service, if the ESSENTIAL PERFORMANCE is lost or degraded due to EM DISTURBANCES.
- PRECAUTION: The performance of the device may be degraded should one or more of the following occur:
 - Operation outside the manufacturer's stated temperature and humidity range.
 - Storage outside the manufacturer's stated temperature and humidity range.
 - Mechanical shock (for example, drop test) or degraded sensor.
 - Patient temperature is below ambient temperature.