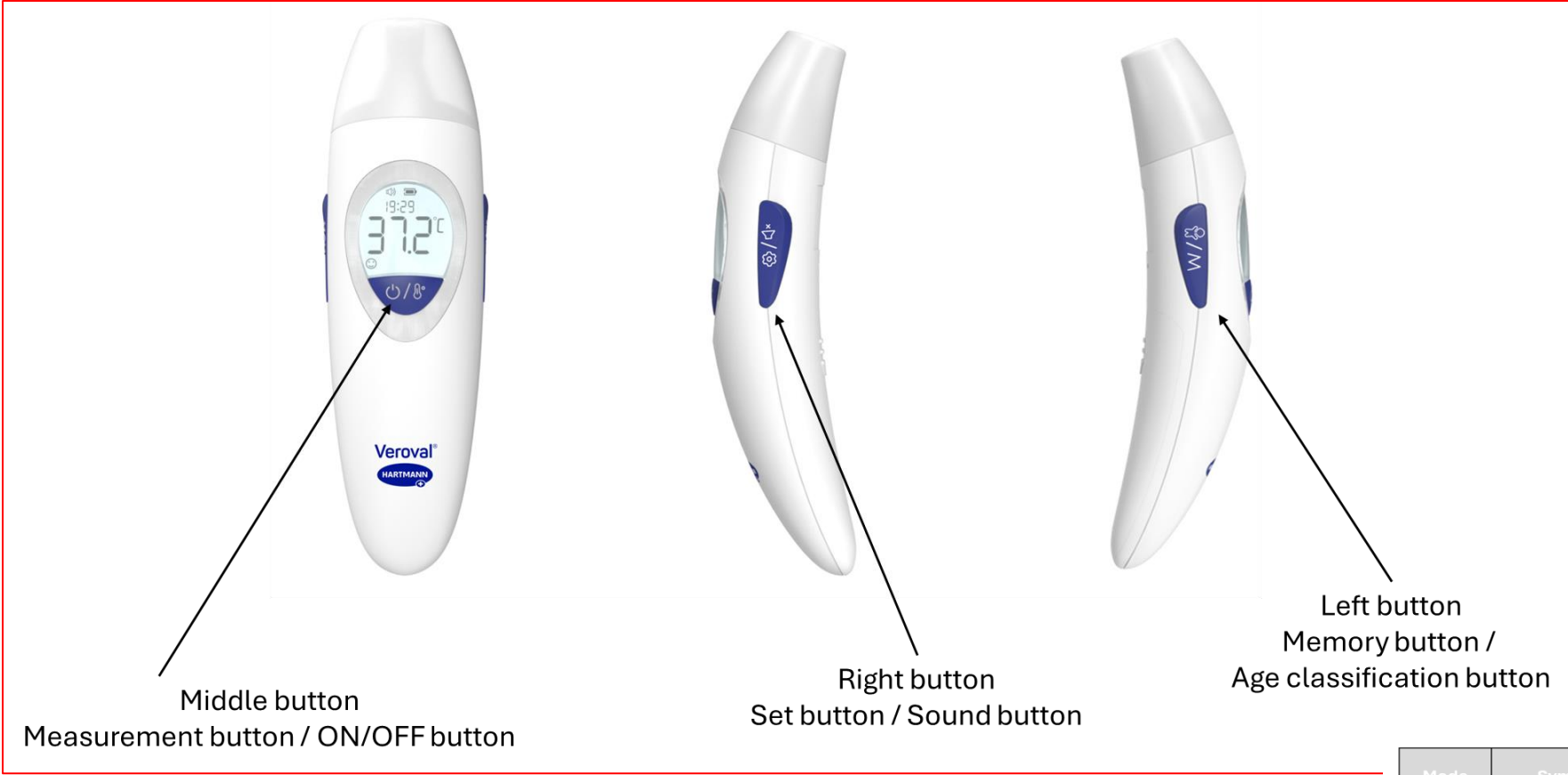




Mode	Symbol
Forehead	
Ear	
Object	

¹ There must be a 5-second interval between two ear measurements. These 5 seconds are needed to allow the reference sensor to cool down and enable a valid measurement.

	Title	Keyflow Hartmann Veroval family – DET-2126			
	Date	2025-01-23	Rev.	03	

- Low battery warning:**
- If the battery load is < 2.70 V, the low battery sign blinks.
 - If the battery load is < 2.60 V, the outer frame of battery symbol is displayed with Lo. With that, no backlight is possible and no measurement can be taken. The device turns off automatically after 30s.

End of measurement:

Object mode – The end of a measurement is signaled through sound/vibration (“-“). The backlight is always green.

Forehead mode/ ear mode – The end of a measurement is signaled through sound/vibration (see table below).

Together with the (alarm) signal/vibration, a colored backlight turns on (green, orange, red). After 15 s, the colored backlight turns off again. For auditive/haptic feedback, the following beep/vibration pattern are implemented:

Temperature [°C]	Temperature [°F]	Beep/Vibration pattern	Backlight color	Smiley indicator	Display
0 – 3 months					
34.0 ≤ T ≤ 37.4	93.2 ≤ T ≤ 99.4	-	Green	😊	
37.4 < T ≤ 43.0	99.4 < T ≤ 109.4	- . . .	Red	😞	
T < 34.0	T < 93.2		Green	😊	Lo
T > 43.0	T > 109.4	- . . .	Red	😞	Hi
3 – 36 months					
34.0 ≤ T ≤ 37.6	93.2 ≤ T ≤ 99.6	-	Green	😊	
37.6 < T ≤ 38.5	99.6 < T ≤ 101.3	-	Orange	😊	
38.5 < T ≤ 43.0	101.3 < T ≤ 109.4	- . . .	Red	😞	
T < 34.0	T < 93.2		Green	😊	Lo
T > 43.0	T > 109.4	- . . .	Red	😞	Hi
36+ months					
34.0 ≤ T ≤ 37.7	93.2 ≤ T ≤ 99.9	-	Green	😊	
37.7 < T ≤ 39.4	99.9 < T ≤ 103.0	-	Orange	😊	
39.4 < T ≤ 43.0	103.0 < T ≤ 109.4	- . . .	Red	😞	
T < 34.0	T < 93.2		Green	😊	Lo
T > 43.0	T > 109.4	- . . .	Red	😞	Hi

Note : “-” long beep (250ms), “.” short beep (180ms)

Automatic switch-off: After 30 s with no button press, the device turns off automatically (without any vibration/sound feedback).

¹ There must be a 5-second interval between two ear measurements. These 5 seconds are needed to allow the reference sensor to cool down and enable a valid measurement.

	Title	Keyflow Hartmann Veroval family – DET-2126		
	Date	2025-01-23	Rev.	03

The key flow described below applies in every "age mode."

Button (left = age/MEM, right = settings, middle = On/OFF)	Button press duration (short: <1s, long: ca. 3s)	Device mode	Sound mode	Measurement mode	Cap mode	Description
☒ left button ☐ middle button ☐ right button	☒ short ☐ long	☐ on ☒ off	☒ on ☒ off (Vibration)	☐ Forehead mode ☐ Object mode ☐ Ear mode	☐ cap on ☐ cap off	Memory mode
☐ left button ☒ middle button ☐ right button	☐ short ☒ long	☐ on ☐ off ☒ on, memory mode	☐ on ☒ off (Vibration)	☐ Forehead mode ☐ Object mode ☐ Ear mode	☐ cap on ☐ cap off	Turn device off by long button press for 4s. Vibration feedback (“-“) signalizes successful turn-off.
☐ left button ☒ middle button ☐ right button	☐ short ☒ long	☐ on ☐ off ☒ on, memory mode	☒ on ☐ off (Vibration)	☐ Forehead mode ☐ Object mode ☐ Ear mode	☐ cap on ☐ cap off	Turn device off by long button press for 4s. Sound feedback (“-“) signalizes successful turn-off.
☒ left button ☐ middle button ☐ right button	☒ short ☐ long	☒ on ☐ off	☐ on ☒ off (Vibration)	☒ Forehead mode ☐ Object mode ☒ Ear mode	☒ cap on ☒ cap off	Device must vibrate if age function mode is changed by button press.
☒ left button ☐ middle button ☐ right button	☒ short ☐ long	☒ on ☐ off	☒ on ☐ off (Vibration)	☒ Forehead mode ☐ Object mode ☒ Ear mode	☒ cap on ☒ cap off	Device must give signal if age function mode is changed by button press.

¹ There must be a 5-second interval between two ear measurements. These 5 seconds are needed to allow the reference sensor to cool down and enable a valid measurement.

	Title	Keyflow Hartmann Veroval family – DET-2126				
	Date	2025-01-23	Rev.	03		
Button (left = age/MEM, right = settings, middle = On/OFF)	Button press duration (short: <1s, long: ca. 3s)	Device mode	Sound mode	Measurement mode	Cap mode	Description
☒ left button ☐ middle button ☐ right button	☒ short ☐ long	☒ on ☐ off	☒ on ☒ off (Vibration)	☒ Forehead mode ☐ Object mode ☒ Ear mode	☒ cap on ☒ cap off	Switching-time between age modes should be the same as in sound-on mode, to ensure consistency between sound-on/sound-off vibration modes.
☐ left button ☒ middle button ☐ right button	☒ short ☐ long	☐ on ☒ off	☐ on ☒ off (Vibration)	☐ Forehead mode ☐ Object mode ☐ Ear mode	☐ cap on ☐ cap off	Turn device on without vibration feedback and when ready to measure “double vibration” as it is in sound mode but here with vibration.
☐ left button ☒ middle button ☐ right button	☒ short ☐ long	☐ on ☒ off	☒ on ☐ off (Vibration)	☐ Forehead mode ☐ Object mode ☐ Ear mode	☐ cap on ☐ cap off	Turn device on without sound feedback and when ready to measure “double beep”.
☐ left button ☒ middle button ☐ right button	☐ short ☒ long	☒ on ☐ off	☐ on ☒ off (Vibration)	☒ Forehead mode ☒ Object mode ☒ Ear mode	☒ cap on ☒ cap off	Turn device off by long button press for 4s. Vibration feedback (“-“) signalizes successful turn-off.
☐ left button ☒ middle button ☐ right button	☐ short ☒ long	☒ on ☐ off	☒ on ☐ off (Vibration)	☒ Forehead mode ☒ Object mode ☒ Ear mode	☒ cap on ☒ cap off	Turn device off by long button press for 4s. Sound feedback (“-“) signalizes successful turn-off.
☐ left button ☒ middle button ☐ right button	☒ short ☐ long	☒ on ☐ off	☒ on ☐ off (Vibration)	☒ Forehead mode ☐ Object mode ☐ Ear mode	☒ cap on ☐ cap off	Swipe-measurement starts (3s), rotating icon on display shows ongoing measurement , sound feedback + 15s backlight when measurement finished, show highest value.

¹ There must be a 5-second interval between two ear measurements. These 5 seconds are needed to allow the reference sensor to cool down and enable a valid measurement.

	Title	Keyflow Hartmann Veroval family – DET-2126				
	Date	2025-01-23	Rev.	03		
Button (left = age/MEM, right = settings, middle = On/OFF)	Button press duration (short: <1s, long: ca. 3s)	Device mode	Sound mode	Measurement mode	Cap mode	Description
☐ left button ☒ middle button ☐ right button	☒ short ☐ long	☒ on ☐ off	☐ on ☒ off (Vibration)	☒ Forehead mode ☐ Object mode ☐ Ear mode	☒ cap on ☐ cap off	Swipe-measurment starts (3s), rotating icon on display shows ongoing measurement , vibration feedback + 15 s backlight when measurement finished, show highest value.
☐ left button ☒ middle button ☐ right button	☒ short ☐ long	☒ on ☐ off	☒ on ☐ off (Vibration)	☐ Forehead mode ☐ Object mode ☒ Ear mode	☐ cap on ☒ cap off	Start measurement, sound feedback + 15 s backlight when measurement finished, after 5 seconds “ready-to-measure” double-beep feedback ¹ .
☐ left button ☒ middle button ☐ right button	☒ short ☐ long	☒ on ☐ off	☐ on ☒ off (Vibration)	☐ Forehead mode ☐ Object mode ☒ Ear mode	☐ cap on ☒ cap off	Start measurement, vibration feedback + 15 s backlight when measurement finished.
☐ left button ☒ middle button ☐ right button	☒ short ☐ long	☒ on ☐ off	☒ on ☐ off (Vibration)	☐ Forehead mode ☒ Object mode ☐ Ear mode	☒ cap on ☐ cap off	Start measurement (1s), sound feedback when measurement finished.
☐ left button ☒ middle button ☐ right button	☒ short ☐ long	☒ on ☐ off	☐ on ☒ off (Vibration)	☐ Forehead mode ☒ Object mode ☐ Ear mode	☒ cap on ☐ cap off	Start measurement (1s), vibration feedback when measurement finished.
☐ left button ☐ middle button ☒ right button	☒ short ☐ long	☒ on ☐ off	☐ on ☐ off (Vibration)	☒ Forehead mode ☒ Object mode ☒ Ear mode	☒ cap on ☒ cap off	Sound on/off, signaled by sound or vibration + icon in display if sound on.

¹ There must be a 5-second interval between two ear measurements. These 5 seconds are needed to allow the reference sensor to cool down and enable a valid measurement.

Title	Keyflow			
	Hartmann Veroval family – DET-2126			
Date	2025-01-23	Rev.	03	

Button (left = age/MEM, right = settings, middle = On/OFF)	Button press duration (short: <1s, long: ca. 3s)	Device mode	Sound mode	Measurement mode	Cap mode	Description
☐ left button ☐ middle button ☒ right button	☐ short ☒ long	☒ on ☐ off	☒ on ☒ off (Vibration)	☒ Forehead mode ☒ Object mode ☐ Ear mode	☒ cap on ☐ cap off	Switch between forehead and object mode.
☐ left button ☐ middle button ☒ right button	☐ short ☒ long	☒ on ☐ off	☒ on ☒ off (Vibration)	☐ Forehead mode ☐ Object mode ☒ Ear mode	☐ cap on ☒ cap off	For ear mode, no action with this kind of button press.
☐ left button ☐ middle button ☒ right button	☐ short ☒ long	☐ on ☒ off	☐ on ☐ off (Vibration)	☐ Forehead mode ☐ Object mode ☐ Ear mode	☒ cap on ☒ cap off	Show °C/°F menu.
☐ left button ☒ middle button ☐ right button	☒ short ☐ long	☐ on ☐ off ☒ off °C/°F-mode	☐ on ☐ off (Vibration)	☐ Forehead mode ☐ Object mode ☐ Ear mode	☒ cap on ☒ cap off	Change between °C/°F.
☐ left button ☐ middle button ☒ right button	☒ short ☐ long	☐ on ☐ off ☒ off °C/°F-mode	☐ on ☒ off (Vibration)	☐ Forehead mode ☐ Object mode ☐ Ear mode	☒ cap on ☒ cap off	Select °C/°F, turn device off (with vibration feedback “-“).

¹ There must be a 5-second interval between two ear measurements. These 5 seconds are needed to allow the reference sensor to cool down and enable a valid measurement.

	Title	Keyflow Hartmann Veroval family – DET-2126				
	Date	2025-01-23	Rev.	03		
Button (left = age/MEM, right = settings, middle = On/OFF)	Button press duration (short: <1s, long: ca. 3s)	Device mode	Sound mode	Measurement mode	Cap mode	Description
☐ left button ☐ middle button ☒ right button	☒ short ☐ long	☐ on ☐ off ☒ off °C/°F-mode	☒ on ☐ off (Vibration)	☐ Forehead mode ☐ Object mode ☐ Ear mode	☒ cap on ☒ cap off	Select °C/°F, turn device off (with sound feedback “-“).
☐ left button ☐ middle button ☐ right button	☐ short ☐ long	☒ on ☐ off	☐ on ☒ off (Vibration)	☐ Forehead mode ☐ Object mode ☐ Ear mode	☒ cap on ☒ cap off	When cap put on/off the device, no vibration feedback.
☐ left button ☐ middle button ☐ right button	☐ short ☐ long	☒ on ☐ off	☒ on ☐ off (Vibration)	☐ Forehead mode ☐ Object mode ☐ Ear mode	☒ cap on ☒ cap off	When cap put on/off the device, no sound feedback.
☐ left button ☒ middle button ☐ right button	☐ short ☒ long (7s)	☒ on ☐ off	☐ on ☐ off (Vibration)	☐ Forehead mode ☐ Object mode ☐ Ear mode	☐ cap on ☐ cap off	Access calibration mode, “.” appears as soon as calibration mode is entered

	Title	Keyflow Hartmann Veroval family – DET-2126		
	Date	2025-01-23	Rev.	03

RELEASE	
Issued / changed by	Checked / approved by
2025-01-23  i.A. Antonia Bernecker Technical Project Manager	2025-01-23  i.A. Anke Reschke Technical Projectmanager

HISTORY	Changes to the previous document will be marked in red and documented below. Before updating the document, all previous red markings have to be set back to black colour.		
Revision	Date	Description	Issued / changed by
01	2024-10-24	Initial version	Antonia Bernecker
02	2024-11-26	Added low battery information; updated end of measurement specification; updated Keyflow according to development decisions	Antonia Bernecker
03	2025-01-23	Updated prolonged sound/vibration duration due to vibration motor specification; added loading-symbol information during forehead measurement; added information on entering calibration mode; updated product picture for button explanation; corrected battery warning voltage	Antonia Bernecker

Beurer Template 410-0005 / rev05 / 2024-01-11